

The Essential  
**LifeStyle** for  
Active, Working, and Traveling  
Oxygen Patients



LIFESTYLE™ PORTABLE  
OXYGEN CONCENTRATOR



Now, for the first time, active, working, and traveling oxygen patients can truly enjoy the lifestyle they desire, with an unrestricted oxygen supply.

## Mobility Redefined

With AirSep's new LifeStyle™ Portable Oxygen Concentrator, those who need the benefits of oxygen therapy can travel and vacation with ease, (including the logistics associated with extended cruises or

airport connections), or simply enjoy unplanned, overnight stays.

Supplemental oxygen users can go wherever adventure calls — without the constant fear of running out of their liquid or cylinder supply. And it's now so easy to pack for a trip to just about anywhere — with one compact, lightweight (less than 10 lb/4.5 kg) unit that makes its own supply of oxygen, whether at home or away.

## No Refills, No Deliveries Ever Again!

The LifeStyle concentrator encapsulates the technology and dependability of a traditional concentrator with the economics of a pulse flow conserving device, providing your patients with five pulse settings, according to prescription. And what's more, you know the peace of mind free from the worry of deliveries and refills, particularly with a highly ambulatory or remotely-located patient.

With LifeStyle, your patient always has an endless supply of oxygen on hand because the device is truly a mini oxygen generating system, producing its own contents everywhere your patient goes.

## Multiple Power Choices

LifeStyle is the most flexible oxygen supply in today's market, powered by normal household AC electric, DC automobile power, or a



rechargeable battery that easily snaps in and out of the unit for quick change-outs.

## Marketing Advantages

Equipment providers and patients who use the new portable oxygen concentrator on a daily basis easily see two primary benefits and a winning combination in LifeStyle: economics and convenience. You can fully promote mobility while eliminating the costly logistic issues of supplying compressed gas or liquid

portable oxygen systems. And patients appreciate all of their oxygen needs met at a fraction of the size and power consumption of traditional concentrators.

The priceless benefits of increased ambulation, essential to a patient's health and well-being, is now within reach.

Endless lifestyle possibilities for the patient and the provider...

## For Active and Traveling Patients...

LifeStyle's multiple power choices enable total flexibility within one versatile, easy-to-use, easy-to-handle unit:

- ▶ AC adapter with universal electrical power supply (100–250 volt, 50/60 Hz)
- ▶ DC adapter (12-volt) for automobiles/motorized boats
- ▶ Rechargeable, internal battery



## ...Or for Leisure Pursuits

Whether going to the beauty parlor, a neighbor's home, a sporting or entertainment event, relaxing with a favorite

*LifeStyle*™  
by AirSep®

Simply change the detachable power cord to LifeStyle's universal AC power supply, as necessary, and you have a unit that can travel with your patients — around the corner or around the globe.



hobby, or just enjoying the convenience of portable oxygen without the fuss of managing contents and rationing supplies, LifeStyle is the answer for today's savvy oxygen user.

## ...For Working Patients

A single, portable, contemporary-styled unit that makes and supplies all of a patient's oxygen needs — all day long — is ideal for use in most full-time or structured work situations.

Even part-time or semi-retired individuals on oxygen therapy can remain gainfully employed or choose to do volunteer activities with a little help from their LifeStyle. LifeStyle is suitable for working at open-air markets, sitting behind the desk in a busy office, cashiering at a store's checkout counter, and a host of other work activities — both indoors and out.



## Convenience at Home...

Even moving about the house has never been easier — with the only oxygen concentrator that patients can actually carry over their shoulder while walking. Climbing up and downstairs while LifeStyle makes its own oxygen — is total confidence for the O<sub>2</sub> patient!



## SPECIFICATIONS

**OXYGEN CONCENTRATION\*:** 1-5 pulse settings; equivalent to a continuous flow of 90% oxygen at  $\pm 3$

**DIMENSIONS:** 5.5 in. H x 7.25 in. W x 16.31 in. L  
(13.9 cm H x 18.4 cm W x 41.4 cm L)

**WEIGHT:** 9.75 lb (4.4 kg)

**POWER:** LifeStyle unit — 13 VDC, 42 Watts  
Power supply — 100-250 VAC,  
50/60 Hz, 2.0 amp maximum  
Battery — 12 VDC, nominal

**BATTERY DURATION:** 50 minutes nominal for fully-charged battery; up to 3-1/4 hours with 3 batteries in optional PowerPack

**BATTERY CHARGE TIME:** 2 to 2-1/2 hours;  
PowerPack with up to 3 batteries —  
4 hours

**AUDIBLE ALARMS AND VISUAL INDICATORS:** Start-up — Audible/visual (green)  
Pulse flow — Visual (green)  
Battery condition — Audible/visual (yellow)  
Cannula disconnect/No inspiration — Audible/visual (red)  
System overdraw — Audible/visual (red)  
High and low pressure — Audible/visual (red)

**WARRANTY:** 2-year parts on LifeStyle unit;  
90 days on rechargeable battery;  
2-year parts on PowerPack

*\*Based on atmospheric pressure of 14.7 psia (101 kPa) at 70° F (21° C)*

## INTERNATIONAL APPROVALS:

LifeStyle Unit –



– CE 0459

– CSA-certified

– AirSep Corporation is ISO 9001:2000-certified




**AirSep Corporation**  
Medical Division

401 Creekside Drive • Buffalo, New York 14228-2085 USA

USA/Canada Toll-free: 800-874-0202

Tel: (716) 691-0202 • Fax: (716) 691-4141

E-mail: [marketing@airsep.com](mailto:marketing@airsep.com) • Internet: [www.airsep.com](http://www.airsep.com)

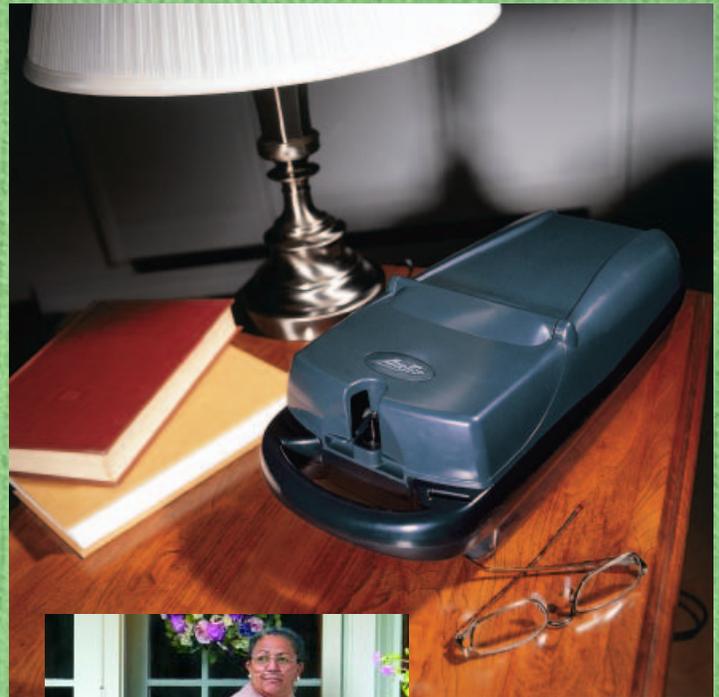

This optional lightweight, retractable cart makes transport of the LifeStyle concentrator and PowerPack (as shown) easy during extended use.



The optional PowerPack is used to recharge LifeStyle batteries and greatly increases ambulation time.

## ORDERING INFORMATION

**PART No. DESCRIPTION**

**AS081-1 LifeStyle Portable Oxygen Concentrator:** includes unit, shoulder strap/carrying bag, rechargeable battery, AC power supply (PW005-1) with detachable power cord, DC automobile adapter (EL204-1), and patient manual.

**Additional travel accessories:**

**BT007-1 Rechargeable battery**  
**BT011-1 PowerPack battery charger**  
**CD023-1 Power cord (120 Volt)**  
**CD023-2 Power cord (120 Volt hospital-grade)**  
**CD017-2 Power cord (220-240 Volt)**  
**CD017-3 Power cord (100 Volt for Japan)**  
**CD025-1 Power cord (220-240 Volt for Australia)**  
**MI210-1 Retractable cart**